Flash&Go
Flash&Go^{Plus}

IMETEC

www.imetec.com

Tenacta Group S.p.A. Via Piemonte, 5/11 24052 AZZANO S. PAOLO (BG) ITALY

Tel. +39.035.688.111 Fax +39.035.320.149

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M1000937
0612 (MMYY)

INTENSE PULSED LIGHT EPILATOR INSTRUCTIONS FOR USE

Dear customer, IMETEC thanks you for purchasing this product. We are certain you will appreciate the quality and reliability of this appliance, designed and manufactured with customer satisfaction in mind. These instructions for use have been drawn up in compliance with the european standard EN 62079.

ATTENTION! Instructions and warnings for a safe use.

Before using this appliance, carefully read the instructions and in particular the safety warnings, which must be complied with. Keep this manual and illustrative guide with the appliance for future consultation. Should you pass the appliance on to another user, make sure to also include this documentation.

Should certain parts of this user manual be difficult to understand or should doubts arise, contact the company before using the product, at the address indicated on the last page.

CONTENTS

| | |
|--|---------|
| Safety warnings | Page 60 |
| Symbols | Page 64 |
| Contraindications for product use | Page 64 |
| Other precautions for safe use | Page 68 |
| Possible side effects | Page 71 |
| Description of the appliance and accessories | Page 74 |
| Treatment areas | Page 75 |
| Hair removal program | Page 76 |
| Information before use | Page 77 |
| Use | Page 81 |
| Maintenance | Page 84 |
| Problems and solutions | Page 86 |
| Disposal | Page 87 |
| Assistance and warranty | Page 87 |
| Skin tone range | I |
| Technical data | II |
| Illustrative guide | III |

SAFETY WARNINGS

- **This appliance cannot be used by persons under the age of 16 and can only be used by persons with reduced physical,**

sensorial or mental abilities or persons lacking experience if they have received prior instructions concerning its safe use and if they are informed about the dangers associated with this product. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

- **If the appliance is used in the bathroom, remove the plug from the socket after use, as bringing it close to water is dangerous even if the appliance is switched off.**
- **To ensure optimal protection, you are advised to install a differential current protective device in the bathroom the nominal differential current of which should not exceed 30 mA (ideally, 10**

mA). Ask the installation technician for advice.

IMPORTANT! Do not use this appliance near bath tubs, wash basins, showers or other water containers! Never place the appliance in water!

- After having removed the appliance from its packaging, verify its integrity based on the drawing and check for any damage caused during transport. If in doubt, do not use the appliance and contact the authorised service centre.
- Do not let children play with the packaging! Keep the plastic wrapping out of the reach of children; risk of suffocation!
- Before connecting the appliance, check that the mains voltage data reported on the identification plate corresponds to that of the available electric mains. The identification plate is attached to the appliance.
- This appliance must be used solely for the purpose for which it was designed, namely, as a intense pulsed light epilator for home use. Any other use is improper and potentially dangerous.
- This appliance is not suitable to be used in the presence of flammable anaesthetic mixtures with air or Oxygen or Nitrous Oxide.
- ⓧ DO NOT use the appliance with wet hands, damp feet or barefoot.
- ⓧ DO NOT pull the power cord or the appliance itself to remove

the plug from the socket.

- ⓧ DO NOT expose the appliance to humidity or weather conditions (rain, sun).
- Always disconnect the appliance from the mains before cleaning or performing maintenance and when the appliance is not in use.
- In case of a fault or malfunctioning of the appliance, switch it off and do not tamper with it. If the appliance needs to be repaired, contact an authorised technical service centre.
- If the power cord is damaged, it must be replaced by an authorised technical service centre, so as to prevent all risks.
- The appliance must only be used on human hair. Do not use on animals.
- ⓧ DO NOT obstruct the ventilation filters (2).
- ⓧ DO NOT use the appliance if the skin colour sensor (5B) is cracked or broken.
- ⓧ DO NOT use the appliance if the external casing is cracked or broken.
- ⓧ DO NOT use the appliance if the lamp (5 or 15) is damaged or if the treatment surface (5A or 15A) is cracked, chipped or missing.
- ⓧ DO NOT use the appliance if the plug (14) is cracked or broken.
- The appliance must only be used with the supplied TYPE KSAPV0361200300D5 plug (14).
- ⓧ NEVER look directly at the light emitted from the treatment surface (5A or 15A).
- ⓧ NEVER point the treatment surface (5A or 15A) whilst emitting pulsed light towards an empty space.

SYMBOLS

Warning.

Double insulation.

Prohibited.

Only for indoor use.

Power supply connector polarity.

CONTRAINDICATIONS FOR PRODUCT USE

Important Safety Information – Read Before Use

This appliance is not intended to be used by just any person. Read and take note of the information in the following paragraph before using the product. For more information and personalised advice visit www.imetec.com or www.bellissima.imetec.com.

If you are unsure whether you can use the appliance or not, ask your doctor or dermatologist!

- **DO NOT USE THE APPLIANCE IF ONE OF THE FOLLOWING CONDITIONS APPLIES:**
- **During pregnancy or breastfeeding.**

- If you have been exposed to the sun or artificial tanning machines in the last 28 days.
- If you have a tattoo or permanent make-up in the treatment area.
- If you have dark brown or black spots, such as large freckles, a birthmark, moles or warts in the treatment area.
- If you have eczema, psoriasis, lesions, open wounds or active infections, such as herpes, in the treatment area. Wait for the affected area to heal before using the appliance.
- If you have a history of keloid formation, if you know you are sensitive to the light (photosensitivity) or take medicine that causes the skin to be more sensitive to light, including non-steroidal anti-inflammatories, (such as aspirin, ibuprofen, paracetamol), tetracyclines, phenothiazines, thiazide, diuretics, sulfonamide, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha-Hydroxy Acids (AHA), Hydroxy Acids (BHA), Retin-A®, Accutane® and/or topical retinoids.
- If you have abnormal skin conditions, such as vitiligo or albinism.
- If you have varicose veins in the areas that are to be treated.
- If you are being treated or have recently been treated with Alpha-Hydroxy Acids (AHA), Beta-Hydroxy Acids (BHA), Retin-A®, topical retinoids or azelaic acid.
- If you have abnormal skin conditions caused by diabetes or other systemic or metabolic disorders.
- If you have been treated with Accutane® (isotretinoin) during the last 6 months.
- If you have been on steroids during the last 3 months.

- If you have a history of herpes outbreaks in the area that is to be treated, unless you have consulted your doctor and received preventive treatment before using the appliance.
- If you suffer from epilepsy.
- If you carry an active system, such as a pacemaker, an incontinence aid, an insulin pump, etc.
- If you suffer from diseases related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- If you have a history of skin cancer or potential skin malignancies.
- If you have undergone radiotherapy or chemotherapy during the last 3 months.
- If you have any other problem for which your doctor does not recommend using the product.

If in doubt, consult your doctor.

THE COLOUR OF YOUR SKIN

The appliance removes unwanted hair and selectively affects the hair pigment. Various amounts of pigment (melanin) are also found in the surrounding skin tissue. The amount of pigment in the skin, which affects the skin tone, determines the risk posed to the skin when any type of pulsed light epilation is used.

Compare your complexion with the skin tone range (A) provided in this manual and/or on the package of the product before using the appliance.

NATURALLY DARK SKIN

DO NOT use the appliance on naturally dark skin!

Using the appliance on dark skin can cause adverse effects such as burns, blisters and changes in the colour of the skin (hyper or hypo-pigmentation).

SKIN COLOUR SENSOR

The appliance has a skin colour sensor (5B) that measures the tone of the skin in the area that is to be treated at the beginning of each session and occasionally during the session. The skin colour sensor (5B) ensures that pulses are only emitted on suitable skin tones.

TANNED SKIN

DO NOT use the appliance on tanned skin or after recent exposure to the sun and/or tanning beds! Such use can cause serious burns, wounds, blisters and changes in the colour of the skin (hyper or hypo-pigmentation).

Tanned skin, especially after sun exposure, contains large amounts of melanin.

This applies to all types of skin and complexions, including those that do not seem to tan quickly. The presence of large amounts of melanin, expose the skin to greater risk when light epilation methods are used.

Avoid exposure to the sun and/or tanning beds for 4 weeks before using the appliance!

THE EYES

NEVER USE the appliance around or near the eyes.

The appliance can be used to remove hair in a delimited part of the face. Special care must be taken to avoid the eyes when treating this area.

Keep away from the eyes!

Do not use the appliance your nipples or genital organs (male or female).

OTHER PRECAUTIONS FOR SAFE USE

BEFORE THE TREATMENT

THE POWER LEVEL

Select the power level CAREFULLY!

The power level refers to the intensity of the pulse of light that is projected on the skin and ranges from the lowest to the highest. The power level LEDs

(12) indicate the set power level. When the power level increases, the results improve, however, the risk of side effects also increases (refer to “Possible side effects”).

Always start the treatment at the lowest power level!

Only if you feel a little or no discomfort with the power level set to the lowest, can you increase the power by one level at the next session, and so on for each subsequent hair removal session. For detailed instructions regarding the power level, refer to “Power Level” in “Information before use”.

- Always shave the area that is to be treated and make sure the skin is clean and dry before using the appliance.
- Cover any birthmark and tattoos before using the appliance.
- Cover dark brown or black spots, such as large freckles, a birthmark, moles or warts before using the appliance.
- ⓧ Never use flammable liquid, such as alcohol, products containing alcohol (including perfume, sanitary products, etc.) or acetone to clean the skin before using the appliance.
- ⓧ Do not use the appliance on any part of the body where you may wish hair to grow eventually.

THE PATCH TEST

It is recommended to perform a patch test before starting treatment, in order to identify the maximum tolerance of your skin to light pulses. Refer to

“Information before use” for instructions regarding the patch test.

DURING TREATMENT

Do not treat the same section of skin more than once per hair removal session!

Do not overlap the pulses!

- STOP USING THE APPLIANCE IMMEDIATELY if blisters form or the skin burns.
- Always make sure that the treatment surface (5A or 15A) is pointed towards the skin and is entirely in contact with it.
- Remove the treatment surface (5A or 15A) from the skin if either (5A or 15A) becomes excessively hot.

AFTER TREATMENT

Do not expose the skin to the sun and/or tanning beds after treatment.

In case of exposure to the sun during the two weeks following a treatment session, protect the treated skin with sun cream SPF 30 or higher.

ATTENTION: treatment may cause a temporary change in pigmentation (refer to “Possible side effects”).

POSSIBLE SIDE EFFECTS

Side effects and complications associated with the appliance are rare when used according to the instructions.

However, any cosmetic procedure, including those intended for home use, pose a potential risk. It is therefore important to understand and acknowledge the risks and complications that may occur when using pulsed light epilation systems designed for home use.

| Negative Event | Probability of the Negative Event | Negative Effect | Probability of the Negative Effect |
|---|-----------------------------------|--------------------------------|------------------------------------|
| Accumulated multiple pulses or overlapping on the same section of skin. | Minimum | Minimal Discomfort to the Skin | Minimum |
| | | Redness on the Skin | Minimum |
| | | Greater Skin Sensitivity | Minimum |
| | | Wounds and burns | Rare |
| | | Scars | Rare |
| | | Changes in Pigmentation | Rare |
| | | Excessive Swelling and Redness | Rare |
| | | Infection | Negligible |
| | | Bruising | Negligible |

Minimal Discomfort to the Skin

Despite home use pulsed light epilation is generally very well tolerated, most users experience mild discomfort during use, which is usually described as a slight feeling of warmth or tingling in the treated sections. This feeling can be felt during the treatment or a few seconds after.

Any other discomfort is not normal and treatment must be stopped because your skin does not tolerate the hair removal treatment or the power level setting is too high.

Redness on the Skin

The skin may become red immediately after using the appliance or within the subsequent 24 hours. The redness usually disappears within 24 hours. Consult your doctor if the redness does not disappear within 2 or 3 days.

Greater Skin Sensitivity

The treated skin becomes more sensitive and may become dry or flaky.

Wounds and Burns on the Skin

Very rarely, burns or wounds may appear after the treatment. Burns or wounds can take several weeks to heal and very rarely, may leave a visible permanent scar.

Scars

Even though it is very rare, permanent scarring may occur. Generally, when scars are formed, they appear flat and white (hypotrophic).

However, the scar may also be big and red (hypertrophic) or big and beyond the wound itself (keloid). Subsequent aesthetic treatment may be required to improve the appearance of the scar.

Changes in Pigmentation

In some rare cases, there lies the risk of temporary hyperpigmentation (increased pigment or brown discolouration) or hypopigmentation (whitening) in the skin surrounding the hair. The risk of changes in skin pigmentation is higher in persons who have a dark skin tone. Generally, discolouration or changes in the skin pigment are temporary, and permanent hyperpigmentation or hypopigmentation is rare.

Excessive Swelling and Redness

In some rare cases, the skin may become very red and swollen. This is more common in sensitive parts of the body. The redness and swelling should lessen in the following 2 to 7 days and ice must be applied often. The area can be cleaned gently, however, avoid exposure to sun.

Infection

Skin infection is very rare, however, it may still occur following a burn or wound caused by the appliance.

Bruising

Very rarely, blue-purple bruising may be caused after using the appliance and may last for 5 to 10 days. As the bruise fades, rust-brown skin discolouration may occur (hyperpigmentation), which may be permanent.

Please consult your doctor or dermatologist in the case of side effects or complications.

DESCRIPTION OF THE APPLIANCE AND ACCESSORIES

Home use pulsed light appliance, based on HPL technology (Home Pulsed Light) to permanently reduce unwanted hair. Permanent hair reduction means long-term reduction in the amount of hairs that grow after a course of treatment.

The appliance emits a light pulse that is absorbed by the pigment in the hair. The light beam is conveyed to the hair follicle, thereby generating heat that can disable it.


Refer to the Illustrative Guide to check the contents of the package. All figures are found in the inside pages of the cover.

1. Control panel.
2. Ventilation filters.
3. Pulse activation button.
4. Power supply connection.
5. Lamp.
 - 5A. Treatment surface.
 - 5B. Skin colour sensor.
6. Ready-for-use LEDs.
 - 6A. Skin colour check LEDs.
7. Error/fault LED.
8. Power level setting buttons (+/-).
9. Stand-by LED.
10. Switch-on/stand-by button .

11. Switch-on LED.
12. Power level LEDs.
13. Lamp replacement LED.
14. Adaptor.
 - 14A. Plug A
 - 14B. Plug B
 - 14C. Plug C
15. Precision Lamp
 - 15A. Treatment surface
 - A. Skin tone range.

For the appliance features, refer to the external packaging. To be used only with the accessories supplied.

TREATMENT AREAS

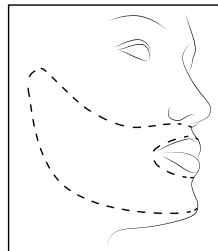The appliance can be used to remove unwanted body hair in the following areas: underarms, bikini line, arms, legs, back and chest.

Women can also use it in the part indicated by the white dashed line.

It is important to only treat the facial areas where there is visible thick and dark superfluous hair. Areas with thin fair hair, often on and under the chin, must not be treated.

HAIR REMOVAL PROGRAM

A typical cycle of hair growth can take 18 to 24 months. During this time, more sessions may be required to achieve permanent hair removal.

Hair removal efficacy varies from person to person and in different body parts, according to the hair colour and level of tolerance of one's own skin.

Recommended treatment plan in a complete hair growth cycle:

- Session 1 to 4: perform the sessions two weeks apart.
- Session 5 to 7: perform the sessions four weeks apart.
- The subsequent maintenance sessions can be performed from time to time, when and if necessary, until long-lasting results are achieved.

Note:

The hair removal results will not be enhanced if the same part of the body is treated more often than recommended (every two weeks).

Dormant hair follicles can be reactivated due to hormonal and physiological changes. Therefore, occasional hair removal sessions may be required.

The individual reaction depends on the type of hair and biological factors that may affect the growth cycle of hair. Some people can react quicker or slower with regards to the recommended hair removal plan.

In any case, a long-lasting result cannot be expected with only one or two sessions.

NOTE

The appliance is not effective on very light hair.

The appliance is not effective on natural white, grey, light blonde and red hair. If your body hair is either of these colours, the appliance will not be effective.

INFORMATION BEFORE USE

POWER LEVEL

The appliance can be set to 5 power levels; 1 (lowest) to 5 (highest), indicated by 3 power level LEDs (12).

SETTING THE POWER LEVEL

Press the power level setting buttons (+/-) (8):

“+” to increase the power level.

“-” to decrease the power level.

The power level LEDs (12) light up as follows according to the set power level:

Level 1 (lowest) - the left LED remains constantly on.

Level 2 - the left LED remains constantly on and the middle LED flashes.

Level 3 - the left LED and the middle LED remain constantly on.

Level 4 - the left LED and middle LED remain constantly on and the right LED flashes.

Level 5 (highest) - all three LEDs remain constantly on.

ATTENTION

It is recommended to always start the first hair removal session by setting the lowest power level (level 1).

If you feel no discomfort during the hair removal session with the power level set to the lowest, you can increase the power to level 2 at the next session, and so on for each subsequent session.

Do not increase the power level if you feel an abnormal discomfort during or after the treatment! (See “Possible side effects”).

Do not increase the power level during a hair removal session, even if no discomfort is felt: side effects can also arise some time after the session.

SKIN COLOUR SENSOR (5B)

Using the pulsed light epilation system on dark skin can cause adverse effects such as burns, blisters and changes in the colour of the skin (hyper or hypopigmentation).

The appliance has a skin colour sensor (5B) that measures the tone of the skin in the area that is to be treated at the beginning of each session and occasionally during the session.

If the skin colour sensor (5B) detects a complexion that is too dark for the treatment, the appliance stops emitting the light pulses automatically and the skin colour check LEDs (6A) start to flash.

If this occurs, try to use the appliance on another part of the body or contact the IMETEC service centre.

COVERING THE ENTIRE TREATMENT AREA

A precise technique must be followed during the session to make sure the entire treatment area is covered and dangerous overlapped pulses are prevented:

- Imagine the treatment area is divided in rows.
- Emit a light pulse proceeding in order, starting from one end of each row and continuing to the other end (from the bottom upwards or right to left, depending on the treatment area).

The area covered by the treatment surface (5A or 15A) leaves temporary marks on the skin due to the light pressure applied by the handpiece.

These marks remain visible for a few moments and are useful to establish the subsequent pulse position.

ATTENTION

Do not overlap the pulses!

Do not treat the same section of skin more than once per hair removal session!
STOP USING THE APPLIANCE IMMEDIATELY if excessive heat or tingling is felt or blisters or small burns appear on the skin!

HOW TO PERFORM THE PATCH TEST

Since the skin tone and sensitivity varies from one body part to another, a patch test should be performed on each part of the body that is to be treated.

- First check the tone of your skin by comparing it to the table on the package and in the instruction manual.

- Identify a small area that is to be treated.
- Shave the area that is to be treated.
- After having set the lowest power level, apply the handpiece to the area and emit the first pulse. If no discomfort is felt, increase the power to the next level, move the handpiece to the adjacent skin area and emit another pulse. Increase the power level on each subsequent area of skin until the maximum tolerable level for the skin type. A feeling of heat can be felt on the skin during the treatment, and this is normal. If discomfort is felt, the maximum tolerable power level is the previous one.
- Check the area after 24 hours: if the redness persists, the maximum tolerable power level for your skin type is the previous one.

Even though it is extremely rare, if hypersensitivity should persist for over hours, it is recommended to consult a doctor.

ATTENTION: the power output directly affects the result, however, due to the principle of heat of the treatment, the maximum limit tolerated by your skin must not be exceeded. This limit is normally high for fair skin and low for dark skin because melanin that darkens the skin increases the production of heat. Do not exceed the maximum limit that your skin tolerates with no discomfort whatsoever.

PREPARING THE AREAS THAT ARE TO BE TREATED

- The skin must be shaved, clean, dry, with no make-up, anti-perspirant or deodorant.
- Moles, birthmarks, dark spots, tattoos and lips must be covered with a white pencil or white patches.


ATTENTION:

Before starting the treatment, read the paragraphs relative to safety warnings, contraindications, precautions for safe use, possible side effects, treatment areas and hair removal program.

USE

PREPARING THE APPLIANCE

1. Remove the handpiece and the other components from the package.
2. Verify that the lamp (5 or 15) is inserted correctly into the appliance.
3. If the appliance is stored in a very cold place, wait for it to reach at least 15°C.
4. If the plug has not been inserted into the adaptor (14), hold the latter, pull the “Open” key downwards and insert the plug (14A, 14B or 14C), which is suitable for the socket of the electrical system, into the relative slot.
5. Insert the adaptor cable (14) into the power supply connection (4).
6. Insert the adaptor (14) in a socket: the stand-by LED lights up (9).

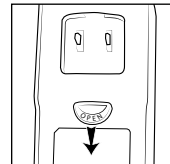TREATMENT

7. Press the switch-on/stand-by button (10): the switch-on LED (11) lights up and the ventilation system is activated (the noise is similar to that of a hair dryer).
8. About 3.5 seconds after having pressed the switch-on/stand-by button (10), the switch-on LED (11) lights up: the appliance is automatically set to the lowest power level.

9. Set the lowest power level when applying the first treatment. Follow the instructions in the “Information before use” chapter if you wish to use a different power level for follow-up sessions.
10. Apply the treatment surface (5A or 15A) firmly on the skin. Only when the treatment surface (5A or 15A) is entirely in contact with the skin, the ready-for-use LEDs (6A) start to flash.
11. Press the pulse activation button (3). The appliance immediately verifies the skin colour: if the skin colour is fair enough for a safe application, the appliance emits the light pulse on the skin together with a sound and the ready-for-use LEDs (6) go off. A slight tingling sensation may be felt. Note: if the appliance does not emit a light pulse and the skin colour check LEDs (6A) flash, the skin tone is too dark for a safe application. Try to use the appliance on another part of the body.
12. Remove the treatment surface (5A or 15A) from the treated area.
13. 3.5 seconds after the previous pulse is emitted, the ready-for-use LEDs (6) light up again.
14. Move the treatment surface (5A or 15A) to another section. Use the pressure marks left previously to place the treatment surface on the adjacent section, avoiding overlapped pulses (See “Covering the entire treatment area”).
15. Repeat the process, starting from point 10, until the entire area to be treated is covered.

AFTER TREATMENT

16. Press the switch-on/stand-by button (10): all the LEDs that are on, go off and the stand-by LED lights up (9).
17. Remember the power level used for the treatment session as this will be reset when the appliance is next switched on.

18. Disconnect the adapter cable (14) from the socket.


Caring for the skin after the hair removal session.

See “Other precautions for safe use”.

TREATMENT WITH THE PRECISION LAMP

If you wish to use the precision lamp (15) for smaller areas or where more precision is required, proceed as follows:

- Before applying the precision lamp (15), emit a light pulse in the area that is to be treated with the standard lamp (5), following the user instructions from 7 to 12. This will allow the skin colour sensor (5B) to verify the tone of the skin that is to be treated and prevent the pulse from being emitted if it is too dark.
- Press the switch-on/stand-by button (10) after emitting a pulse: the switch-on LED (11) goes off and the standby LED (9) goes on.
- Remove the lamp (5) from the handpiece by grasping and pulling it gently upwards to release it.
- Insert the precision lamp (15) into the handpiece by pressing gently.
- Follow the user instructions from point 7 to point 18, being careful not to go over the area previously covered by the emitted pulse of the standard lamp (5).
- The procedure must be repeated every time the power cord (14) is disconnected from the socket.

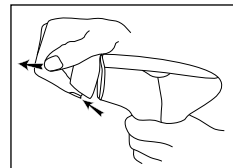MAINTENANCE

CLEANING AND STORING

It is recommended to clean the appliance and the treatment surface (5A or 15A) after each hair removal session:

- Disconnect the appliance from the adaptor (14) before cleaning it.
- Clean the surface of the handpiece and the treatment surface (5A or 15A) gently, with a dry, clean cloth.
- Store the appliance in a cool dry place, away from heat and water.


ATTENTION:

Never place the appliance or any of its parts in water!

REPLACING THE LAMP (5 or 15)


- When 100 available light pulses remain, the the lamp replacement LED (13) starts to flash.
- When the pulses are about to end, the lamp replacement LED (13) starts flashing faster.
- When no more lamp (5 or 15) pulses are left, the lamp replacement LED (13) remains on and no further pulses can be emitted.
- The lamp (5 or 15) must be replaced with an original one for the relative model in order to continue the treatment session.

Note: do not use the lamp and have it replaced if large stains appear inside the lamp (5 or 15) or the treatment surface (5A or 15A) is cracked or broken.

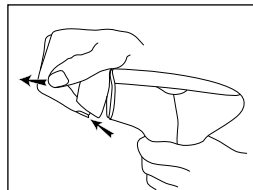1. Press the switch-on/stand-by button (10) to set the appliance in Stand-by mode.
2. Disconnect the adaptor (14) from the socket.
3. Remove the lamp (5 or 15) from the handpiece by grasping and pulling it gently upwards to release it.

4. Unwrap a new lamp (5 or 15).
5. Insert the new lamp (5 or 15) into the handpiece by pressing gently.

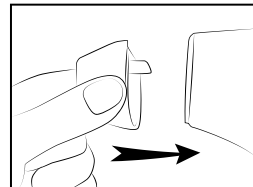Note

After having replaced the lamp (5 or 15) with a new one, always set the power level to a lower one than that used normally. Increase the power level for the subsequent session only if no abnormal discomfort is felt during or after the treatment.

PROBLEMS AND SOLUTIONS

| | |
|--|---|
| <p>The appliance does not go on</p> | <ol style="list-style-type: none"> 1. Make sure the adaptor (14) is inserted correctly into the power supply connection (4). 2. Make sure the adaptor (14) is inserted into a wall socket. |
| <p>No light is emitted when the pulse activation button (3) is pressed.</p> | <ol style="list-style-type: none"> 1. Make sure the ready-for-use LEDs (6) flash, which indicates good contact with the skin and the surface treatment (5A or 15A) is evenly and firmly pressed on the skin. OR 2. Check the skin colour check LEDs (6A): if they flash, the skin tone is too dark for a safe application. Try to use the appliance on another part of the body or contact the IMETEC service centre. OR 3. Check the lamp replacement LED (13): if it is on, replace the lamp (5 or 15) by following the instructions in "Replacing the lamp". OR 4. Verify that the ready-for-use LEDs (6) are on. If the ready-for-use LEDs (6) remain off after 10 seconds, switch the appliance off and on again by pressing the switch-on/stand-by button (10) twice. If the problem persists, contact the IMETEC Customer Service Centre. OR 5. Check the error/fault LED (7): if the LED is on, switch the appliance off and on again by pressing the switch-on/stand-by button (10) twice. If the problem persists, contact the IMETEC Customer Service Centre. |

DISPOSAL

The packaging is made up of recyclable material. Dispose of it in compliance with the environmental protection standards.

The appliance must be disposed of properly, in accordance with the provisions of European standard 2002/96/EC.

The recyclable materials contained in the appliance must be recovered so as to prevent environmental pollution.

For further information, contact the local disposal body or the appliance dealer.

ASSISTANCE AND WARRANTY

To request repairs or purchase spare parts, please contact IMETEC's authorised customer service centre on the Toll-Free Number below or via the website. The appliance is covered by the manufacturer's warranty.

For details, read the attached warranty sheet. Failure to comply with the instructions for use and maintenance provided in this manual will render the manufacturer's warranty null and void.

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